
ELECTRIC BEAUTY CHAIR

PRESTIGE GOLD

MANUAL

MANY BEAUTY

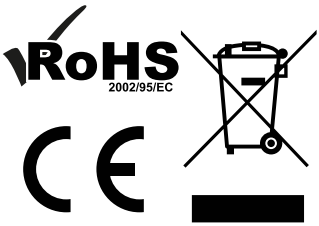


Responsible entity:
MANYBEAUTY Radosław Juszczak
ul. Młyńska 61/2
75-422 Koszalin, Polska
NIP: 669-239-72-28
E-mail: sklep@manybeauty.pl
www.manybeauty.pl

EU Compliance: EMC, LVD,
Machinery Directive
2006/42/WE

1. General Information
Manufacturer and responsible entity:
MANYBEAUTY Radosław Juszczak
ul. Młyńska 61/2
75-422 Koszalin, Poland
VAT ID: 669-239-72-28
E-mail: sklep@manybeauty.pl
Website: www.manybeauty.pl

Model: Prestige Gold
Intended use: For professional use in beauty salons, aesthetic medicine clinics, and SPA centers.
Maximum load capacity: 250 kg
EU Compliance: CE, RoHS, EMC, LVD, Machinery Directive 2006/42/EC
Country of manufacture: China

2. Set Contents
- Chair with electric base and PVC upholstery
- Removable armrests x 2
- Leg rest + extendable foot support
- Headrest with face hole + removable cushion
- Velvet back cushion
- Hand-held remote control
- Built-in dual foot control
- Adjustable leveling feet x 6 + 2 spares
- Power cable (108 cm, EU plug)
- User manual

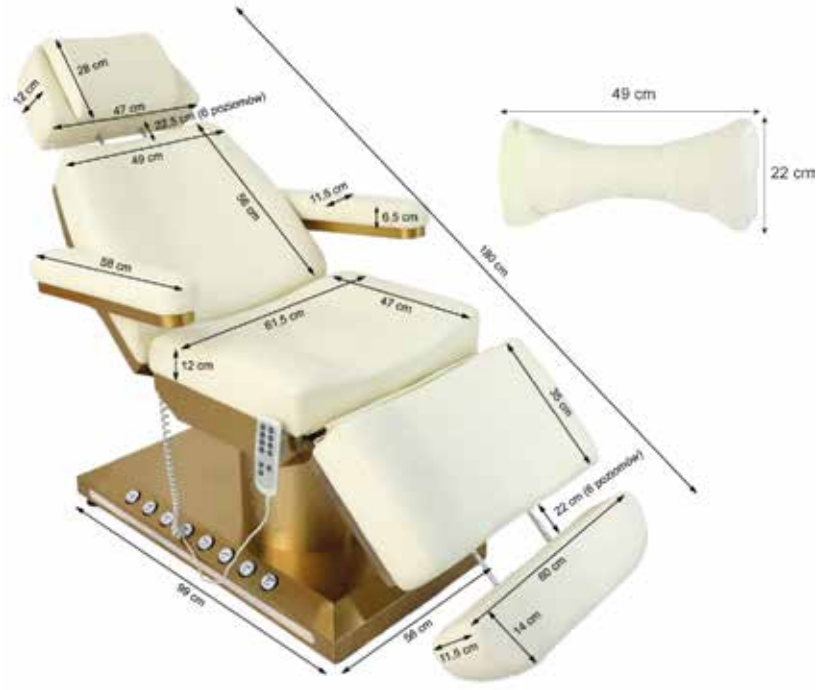

Parameter	Value
Power	200 W
Net weight	80 kg
Package dimensions	131 × 63 × 66 cm
Upholstery color	Ecru
Base color	Gold
Height adjustment range	66–90 cm
Total length	180 cm (over 200 cm when fully extended)
Seat width	61.5 cm
Width with armrests	84.5 cm
Base dimensions	56 × 99 cm
Foam thickness	12 cm
Motors	4 electric motors

4. Functions and Adjustments

- Height adjustment: 66 – 90 cm
- Backrest adjustment: 0 – 90°
- Leg rest adjustment: 0 – 78°
- Seat tilt adjustment
- LED lighting in the base
- Hand-held remote and dual foot control

5. Assembly and First Use

- Carefully unpack the chair, paying attention to delicate parts and finish.
- Install the leveling feet (6 pcs + 2 spares) on the bottom of the base.
- Connect the power cable to a socket compliant with EU standards (230 V ~ 50 Hz).
- Attach the armrests, leg rest, and headrest by inserting them into the designated sockets.
- Test the control system using the hand-held remote and foot switches.
- Test each function separately before first use with a client.

6. Operating Instructions

- Before each use, check the chair's stability.
- Ensure the client is seated or lying evenly – avoid one-sided overloading.
- Do not exceed the maximum load of 250 kg.
- Do not sit or stand on moving parts (e.g., the leg rest).
- Use only indoors, in dry environments.
- Wipe surfaces with a damp cloth and a PVC-safe disinfectant.

7. Maintenance and Cleaning

- Clean upholstery with a soft cloth and a mild detergent suitable for PVC surfaces.
- Do not use corrosive agents or alcohol.
- Regularly check the functionality of the remote and cables.
- If any structural looseness is observed – immediately contact the responsible entity.

8. Warnings and Safety Precautions

- Electrical device – do not use near water!
- Do not open the motor or control panel casings!
- All warranty repairs must be carried out by the responsible entity or an authorized service provider.
- Keep the device away from children.
- Disconnect the chair from power after use.
- Before transport or moving, disconnect from power and secure all movable parts.

9. Disposal

This product is marked with the symbol for selective collection of electrical and electronic waste. Do not dispose of it with regular household waste. Deliver it to an appropriate electronic waste collection point in accordance with the WEEE Directive (2012/19/EU).

10. Warranty and Service

The product is covered by a 12-month manufacturer’s warranty.

All complaints should be sent to: sklep@manybeauty.pl with a description of the issue and proof of purchase.

For many other products visit

www.manybeauty.pl